

Food Screen Assay Card

This reagent kit is designed to test for *Campylobacter*, *Listeria monocytogenes*, and *Salmonella*, from two unknown liquid or dry samples. These protocols are to be performed manually and have been simplified to reduce the risk of operator error.

WARNING: If the unknown sample to be tested is suspected to contain pathogen or toxin, it must be handled by trained personnel and treated as if it were highly poisonous or infectious throughout the sample purification procedure, the assay procedure, and when disposing of all waste materials that have come in contact with the suspected sample.

Hybridization Probe Assays

Kit Contents

1- 5 mL Reagent Grade Water	2- 5 mL Unknown Sample Bottles	4- 1 mL Syringes with Cannula Tips
1- Sample Swab Pack	1- Instruction Booklet	1- Freeze-dried Reagents in a Foil Bag
1- Transfer Pipettes Pack	1- Loading Instructions Sticker	

Kit Part No. PATH-ASY-0010

Part # RAZR-ASY-0003, RAZR-ASY-0013, and/or RAZR-RED-4000 may be required

Programming the Run Protocol

Before the sample is loaded into the pouch, verify that the run protocol for this test is loaded onto the instrument. If it is not and you are using a RAZOR®, use the RAZOR desktop software to create the protocol and load it onto the instrument. Use the data in this card when creating the protocol for this particular test with the Protocol Wizard. For the RAZOR EX, you can load the protocol by scanning the square protocol bar code on the reagent box. If the bar code is damaged, use the generic bar codes that are located below. Times, temperatures, and results analyses are all included in the protocol bar code.

Loading Protocols for the RAZOR® EX

After the square protocol bar code has been scanned, scan the rectangular bar code on the pouch fitment. If the rectangular bar code is damaged, use the generic one below. **Note: you will only be able to use the generic bar code for one run. In order to re-use it, you have to delete the run.**



FOODSCNA-00000001



Protocol Code: FOODSCNA

Step 1:

Negative Port

Add 0.5 mL reagent grade water

Step 2:

Unknown 1 Port

Add 0.5 mL sample 1

Step 3:

Unknown 2 Port


Add 0.5 mL sample 2

Step 4:

Positive Port

Add 0.5 mL reagent grade water

Rectangular Bar Code

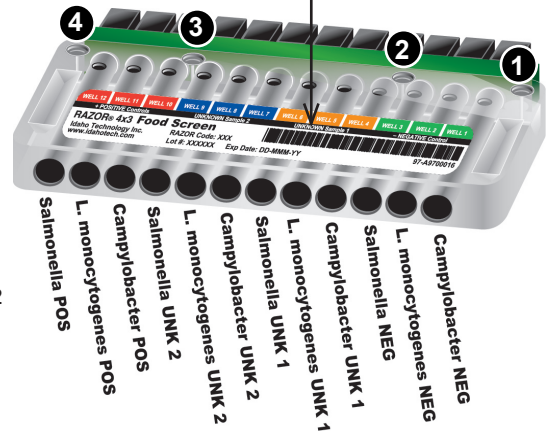

Protocol Steps for the RAZOR®

Protocol Steps	Data to be Entered						
Step 1: Protocol Identifier	Enter: 0AH to identify this protocol. This identifying code matches the first three characters on the pouch. This will help you identify the protocol after it is loaded on the instrument.						
Step 2: Protocol Parameters	Initial Hold Temperature: 93°C	Denature Temperature: 91°C	Anneal Temperature: 60°C				
	Initial Hold Time: 120 s	Denature Hold Time: 3 s	Anneal Hold Time: 15 s				
	Cycles: 55						
Step 3: Select Organisms for the Protocol	A list of organisms is listed in the All Organisms column. Select Campy and click the Add button. Repeat for L. mono and Salm . All three have just been added to the Organism For This Protocol column. If Campy , L. mono or Salm is not listed, see the <i>RAZOR Pouch Instruction Booklet</i> for instructions on how to add new organisms.						
Step 4: Add Environmental Samples	Add two unknown samples by clicking the Add button and entering UNK1 . Repeat the steps and enter UNK2 .						
Step 5: Assign an Organism to Each Sample Position	Assign Campy to positions 1, 4, 7 and 10 Assign L. mono to positions 2, 5, 8 and 11 Assign Salm to positions 3, 6, 9, and 12						
Step 6: RAZOR Pouch Definition	Fill out the definition as shown in the image.	Razor Bag Definition					
		Pos	Organism	Target	Type	Env. Sample	Concentration
		1	Campy	Target	NEG		
		2	L. mono	Target	NEG		
		3	Salm	Target	NEG		
		4	Campy	Target	UNK	UNK1	
		5	L. mono	Target	UNK	UNK1	
		6	Salm	Target	UNK	UNK1	
		7	Campy	Target	UNK	UNK2	
		8	L. mono	Target	UNK	UNK2	
		9	Salm	Target	UNK	UNK2	
		10	Campy	Target	POS		
		11	L. mono	Target	POS		
		12	Salm	Target	POS		
Step 7: RAZOR Meta Calls	This step displays an image of the conditions required for an unknown to be called positive. You do not need to enter any data in this step. Click Finish to save the protocol.						
NOTE:	Before interpreting results on machine's LED screen, remember assays correlate to the channels shown in diagram to right.						

